


DIOLASHeerIce - DL1200 mini

Latest Diode Laser Hair Removal



- Comfortable
- Safety
- Effective
- Fast



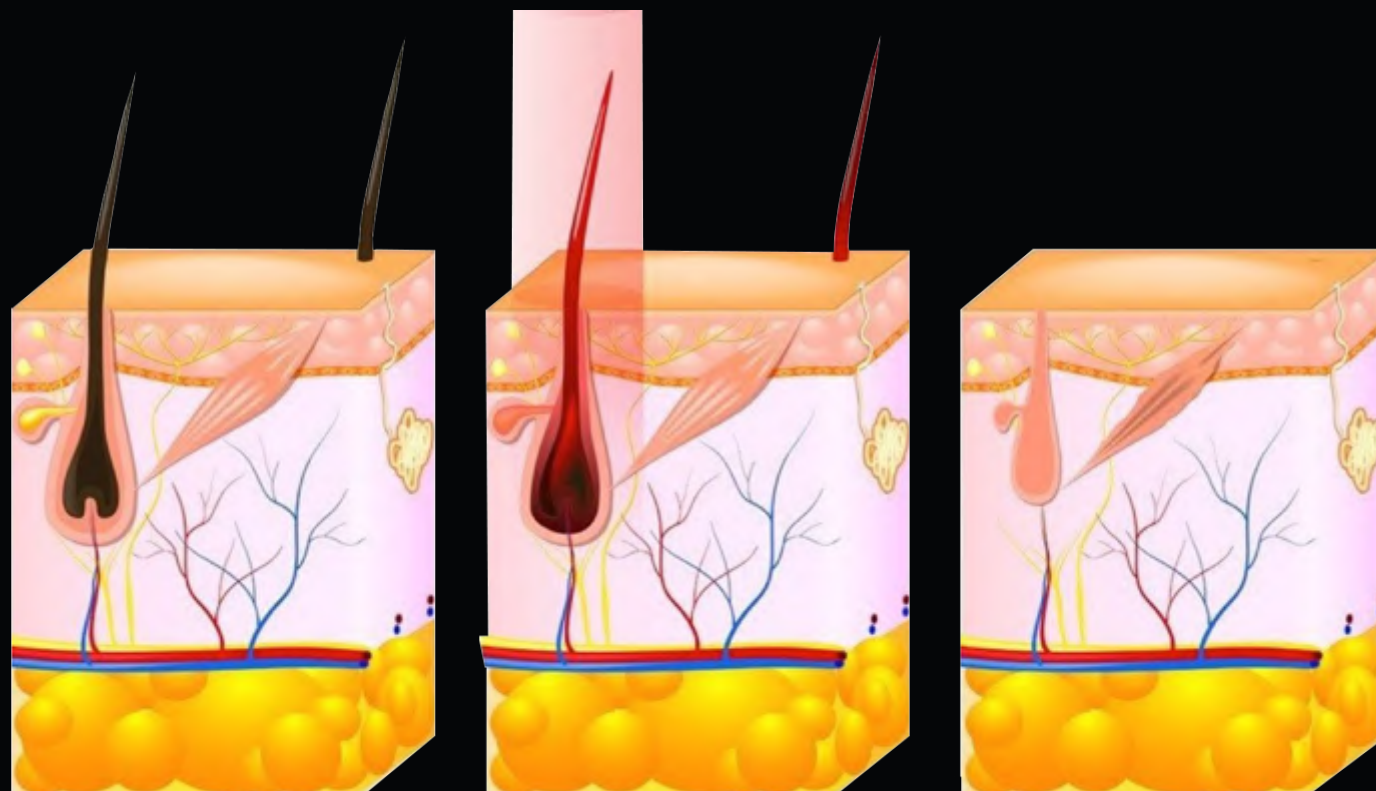


Before

The melanin absorbs
the laser and heats up

After

PRINCIPLE OF TREATMENT

The principle of diode laser hair removal treatment is based on selective photothermal decomposition. The melanosomes in the hair follicle can selectively absorb the energy of the laser. The laser energy emitted by the machine is easily absorbed by the colored hair follicle without damaging the epidermal tissue. The energy emitted by the laser is absorbed by the hair and the pigment in the hair follicle. It is converted into heat, thus increasing the temperature of the hair follicle. When the temperature rises to a certain level, the hair follicle will be irreversibly damaged, and the hair will lose its original environment and be completely removed.

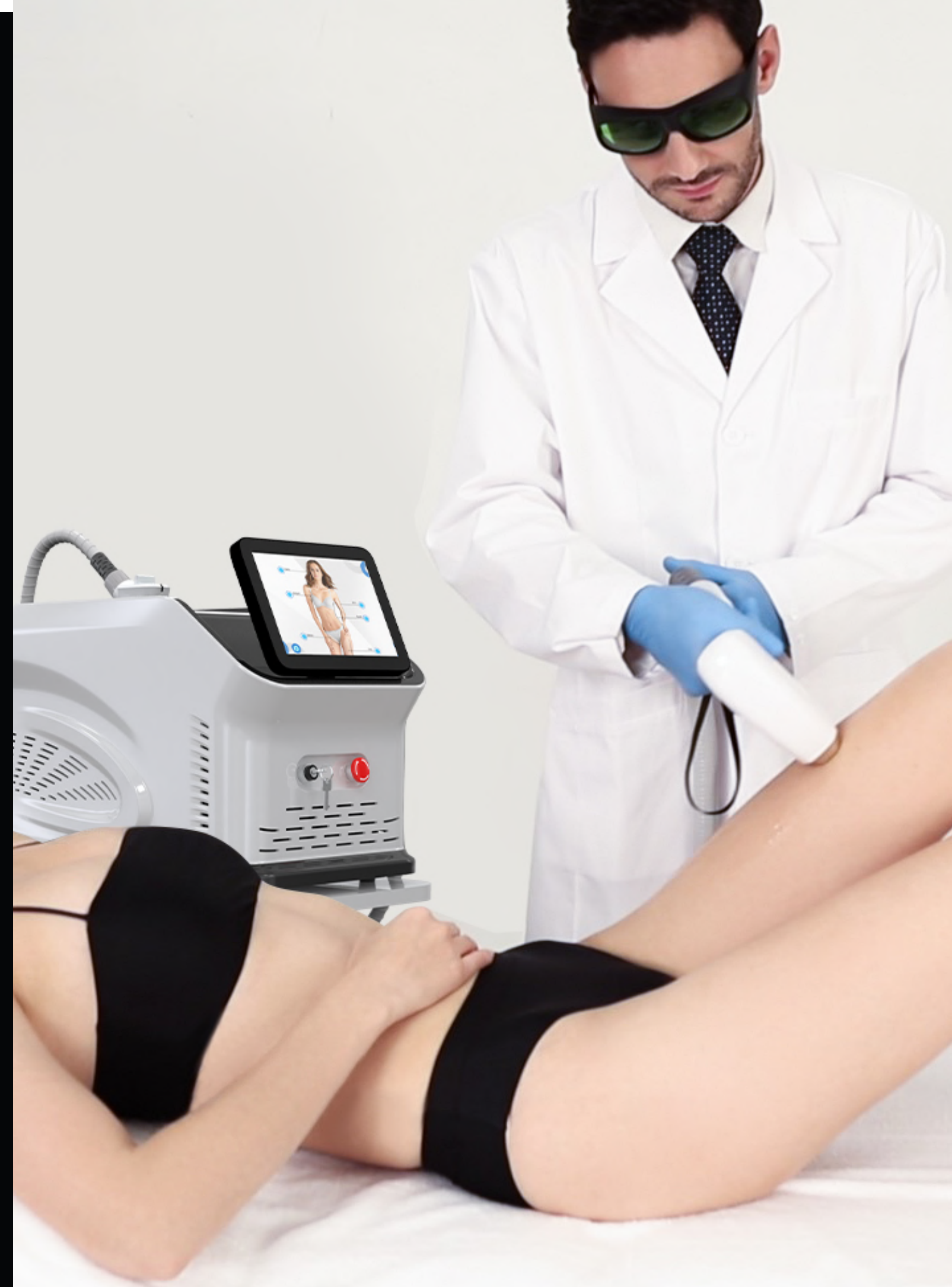

TRIPLE WAVELENGTHS, BEST RESULT FOR ALL SKIN TYPES AND ALL HAIR TYPES




755nm offers more powerful energy absorption by the melanin chromophore, making it ideal for light colored and thin hair. It targets the Bulge of the hair follicle and is especially effective for superficially embedded.



808nm is classic wavelength in laser hair removal, the 808nm wavelength, offers deep penetration of the hair follicle with high average power, suitable for most skin types and hair types.




1064nm wavelength is characterized by lower melanin absorption, making it a focused solution for darker skin types. At the same time, the 1064nm offers the deepest penetration of the hair follicle.

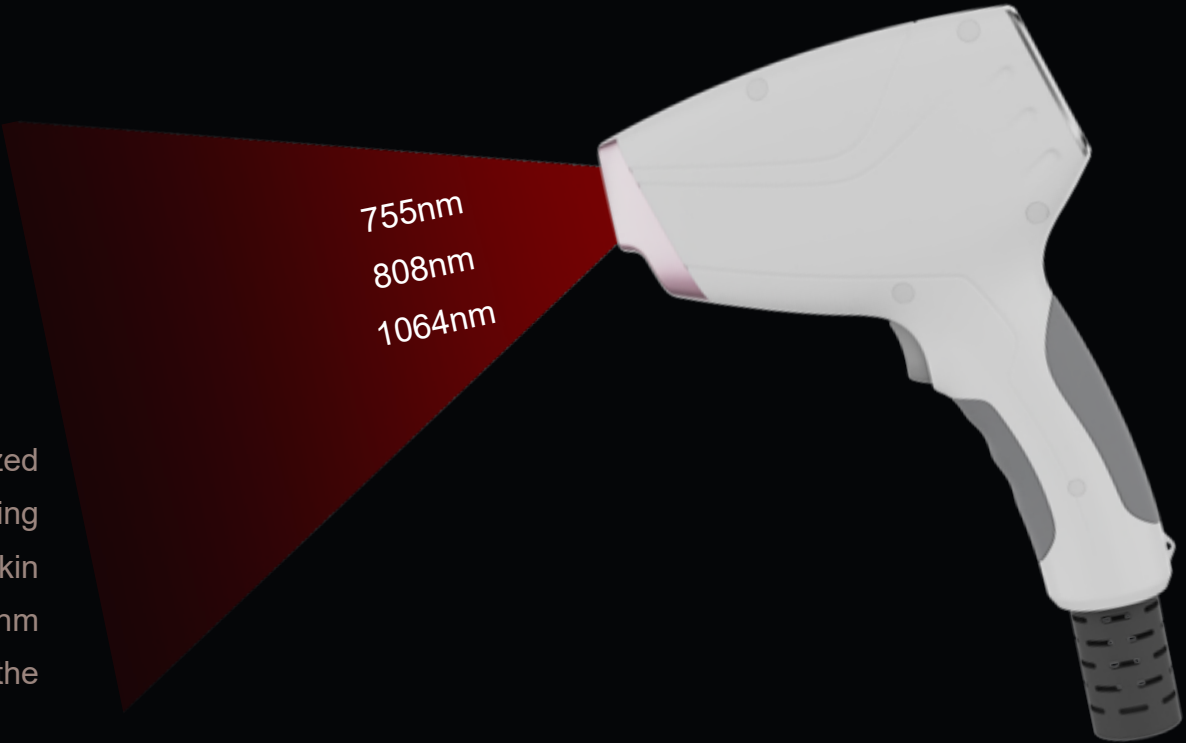




600W-2000W MULTIPLE POWER COMBINATIONS OPTIONAL


A variety of different powers and different spot sizes are available to meet the treatment needs of different parts of the customer

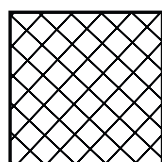


12*12mm
600W 1200W

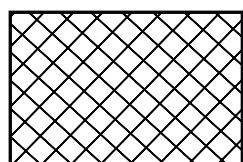

12*20mm
1200W

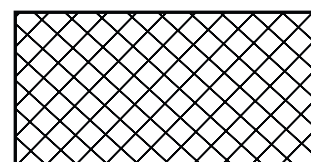

12*24mm
1600W 2000W

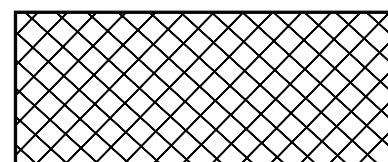

12*28mm
2000W

FACIAL TIP

Reaching the hard to reach


The special Facial TIP allows treatment of hard to reach areas, including the ears, nostrils and glabella.



20HZ ULTRA-FAST FREQUENCY

The maximum frequency can reach 20 Hz, which means that laser is emitted 20 times per second. The treatment speed is twice as fast as other companies' products.

AMERICAN COHERENT LASER BAR

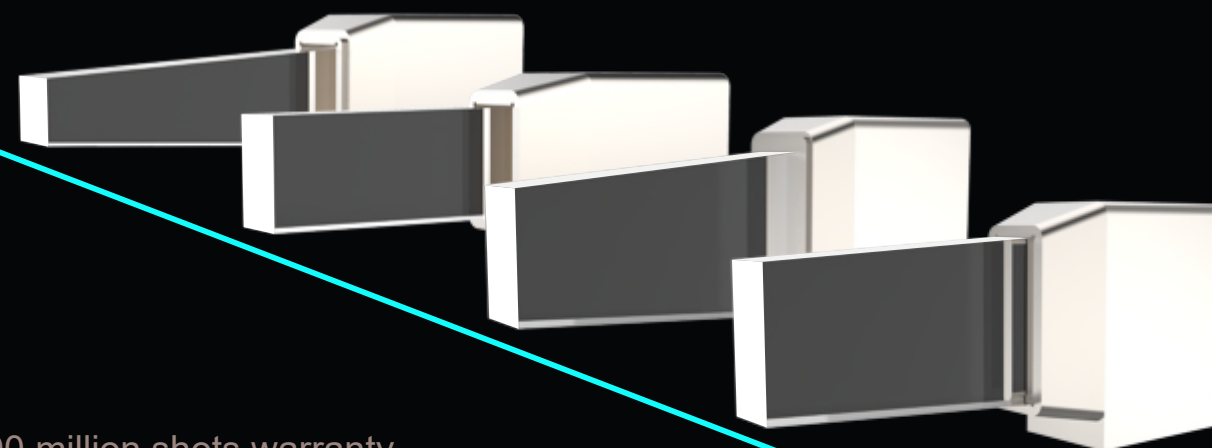

100 million shots warranty

Crystal and laser are packaged in one, with low energy loss and long life

TOUCH SCREEN HANDLE

More than just display

1.54-inch high-definition OLED touch screen, treatment parameters can be adjusted directly on the handle screen.



Energy adjustment



Pulse width adjustment




Frequency adjustment



Real-time temperature monitoring



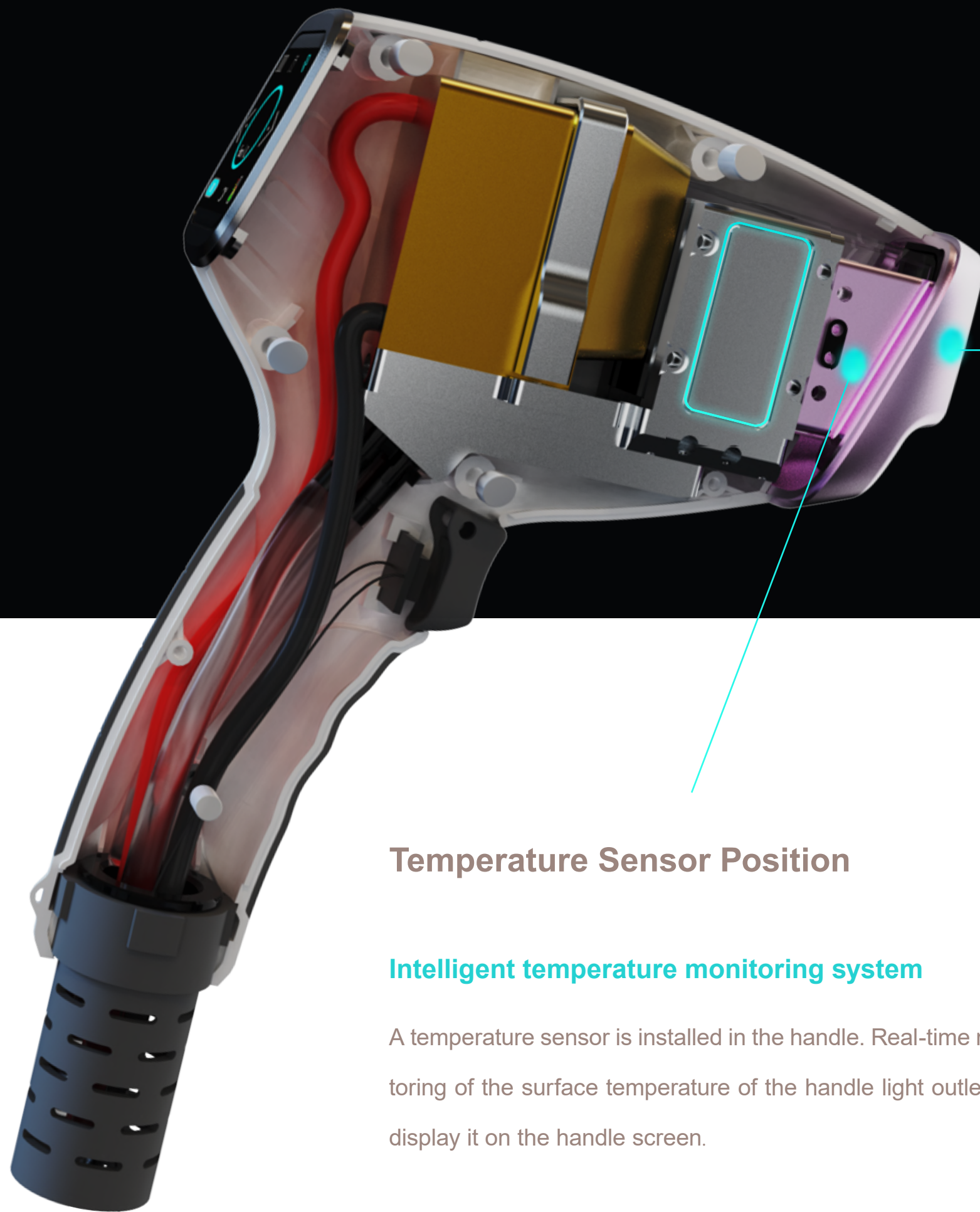

Sapphire crystal

Best light and cold conductivity, ensuring the effectiveness and comfort of the treatment



Temperature Sensor Position

Intelligent temperature monitoring system

A temperature sensor is installed in the handle. Real-time monitoring of the surface temperature of the handle light outlet and display it on the handle screen.

Super refrigeration system

Two semiconductor cooling plates, 108w super cooling power, the surface temperature of the spot can reach -30 °C



INTELLIGENT OPERATING SYSTEM

The operator only needs to select the part and skin color to be treated, and the system will automatically generate the recommended treatment parameters. Even an inexperienced operator can easily grasp the operation of the instrument.

7 Different Languages Interface Available



20 REAL-TIME INTELLIGENT MONITORING

The whole machine is equipped with 20 intelligent detection functions, including water flow detection, water temperature detection, handle temperature detection, handle connection status detection, water level detection, cooling sheet work detection, etc.



Handle connection status detection

Cooling sheet work detection, etc


Water level detection

Water flow detection

Handle temperature detection

Particulates
0.01-0.1μm

Metal Ion

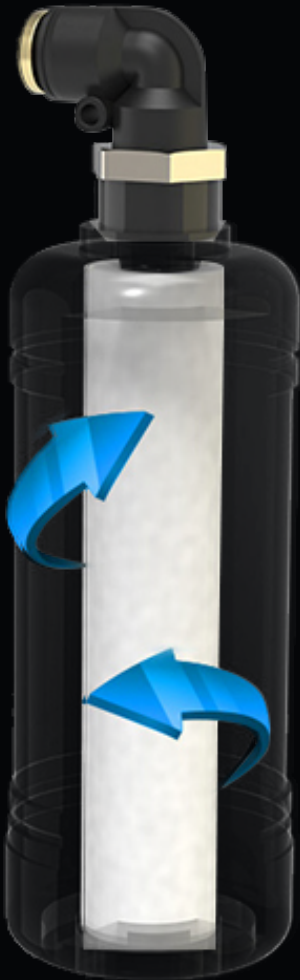

PP filter



Ion filter

2 FILTERS, QUADRUPLE PROTECTION

Each handle is fitted with an advanced two-stage water filtration system. The first stage adopts PP cotton to filter out impurities and prevent laser blockages, the second stage uses a special Ion filter to filter out metal ions, avoiding inner laser corrosion and prolonging system life.

DIOLASHeer ice



SPECIFICATIONS

Power	3000W
Power of handle	600-2000W optional
wavelength	808nm (755nm+808nm+1064nm optional)
Machine screen	12.1 inch
Handle screen	1.54 inch
Energy density	1-120J/cm2 (Deviation $\leq\pm 2\%$)
Pulse width range	1-200ms
Spot size	12*12mm ; 12*20mm; 12*24mm; 12*28mm
Frequency	1-10HZ (600-1200w) 1-20HZ (1600-2000w)
Cooling System	TEC cooling system
Net Weight	40kg
Dimension	400*530*440mm
Package Size	740*580*700mm
Fuse specification	Ø5×25 10A
Voltage	AC220V $\pm 10\%$ 10A 50HZ , 110v $\pm 10\%$ 10A 60HZ