

SMAS FOCUS V

/ Skin Tightening · Face Lifting · Body Slimming · Vaginal Tightening /


5D HIFU



One Machine, Multi-Functional

Only need to replace the handle, you can upgrade the product function at any time.

How it works

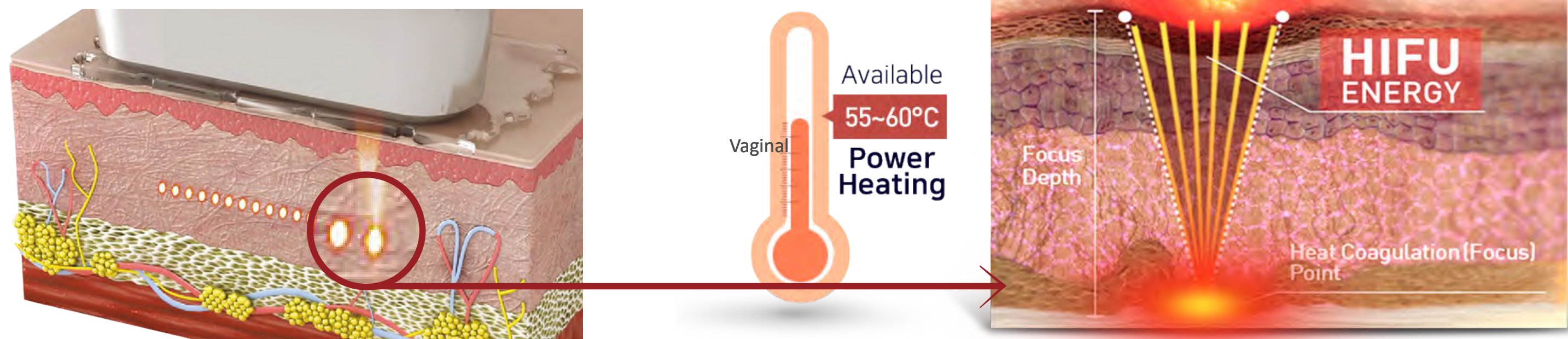

HIFU-powered transducers deliver concentrated waves of ultrasound into multiple layers of skin and body tissues without making direct contact with surrounding areas.

Transferring stabilized temperatures of 65~75°C into the targeted zones, the SMASFOCUS V aims to induce a process of neocollagenesis in the face while tightening body tissues for contouring purposes


FOR LIFTING & TIGHTENING



1 The dermis layer of the skin naturally synthesizes less collagen over the course of time.



3 Stimulated by the focused energy, skin collagen starts to regeneration.

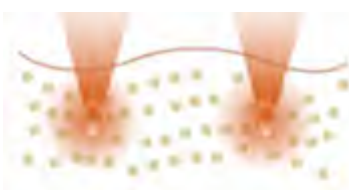


2 Focused ultrasound waves deliver energy directly within the skin layer causing a shrinking effect.



4 Through the regeneration of collagen, wrinkles are eliminated in the process of both firming and tightening the skin.


FOR SLIMMING & CONTOURING

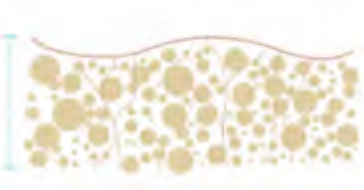

1 Fat cells and cellulite is formed within the body.



3 As a natural healing response of the body the cells regenerate joining together to fill the gap.

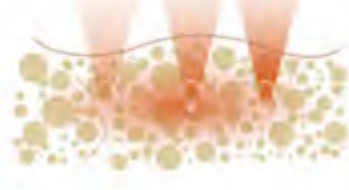


2 The SMASFOCUS V energy destroys the fat cells that have been targeted.



4 In result, the layer of fat naturally exits the body delivering a smooth, contoured body line.

FOR VAGINAL TIGHTENING

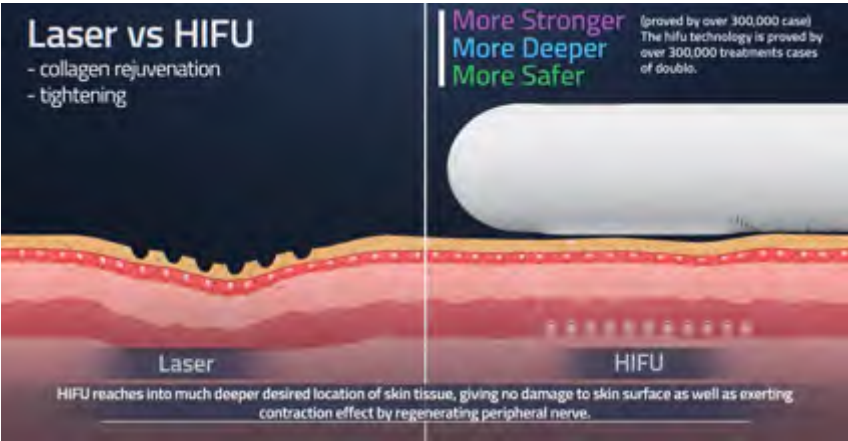

Product Introduction





Vmax single dot handle



1.5mm, 3.0mm, 4.5mm

Eye brow, Wrinkles around eyes,
Wrinkles around month



2D hifu handle



1.5mm, 3.0mm, 4.5mm

Check, Forehead, Double chin,
Neck wrinkles, V-line forming



3D hifu handle



1.5mm, 3.0mm, 4.5mm

Check, Forehead, Double chin,
Neck wrinkles, V-line forming



6.0mm, 8.0mm, 10.0mm, 13.0mm, 16.0mm

Abdomen, Inner thighs, Outer
thighs, Hips & Buttocks, Inner
arms, Breast



Liposonix handle



8.0mm, 13mm

Body slimming and fat reduction



Vaginal handle



3.0mm, 4.5mm

Vaginal tightening

Forehead Wrinkle

Eye brow

Crow's Feet

Cheek

Skin Tone
Improvement

Nasolabial Folds

V-line Forming

Jowl Line

Double Chin

Neck Wrinkles



Easy to upgrade


Each machine is pre-loaded with a software system with all handle functions before shipment. When users need to upgrade the machine functions, they only need to purchase the handle with that function, plug the handle into the machine, and the software system can automatically identify the new handle to realize machine function upgrade.




2D HIFU



3D HIFU

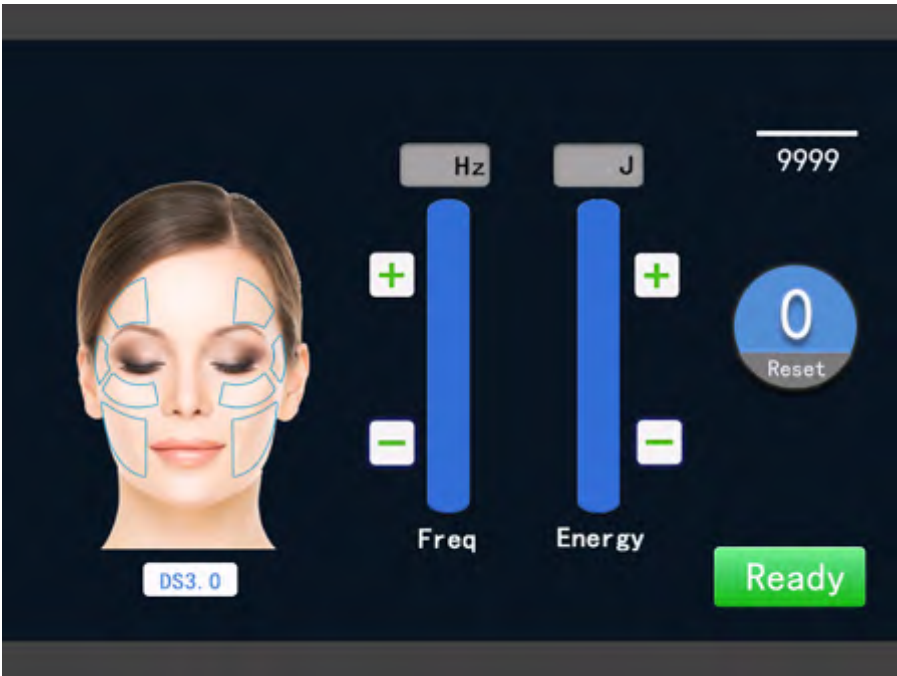

Vmax

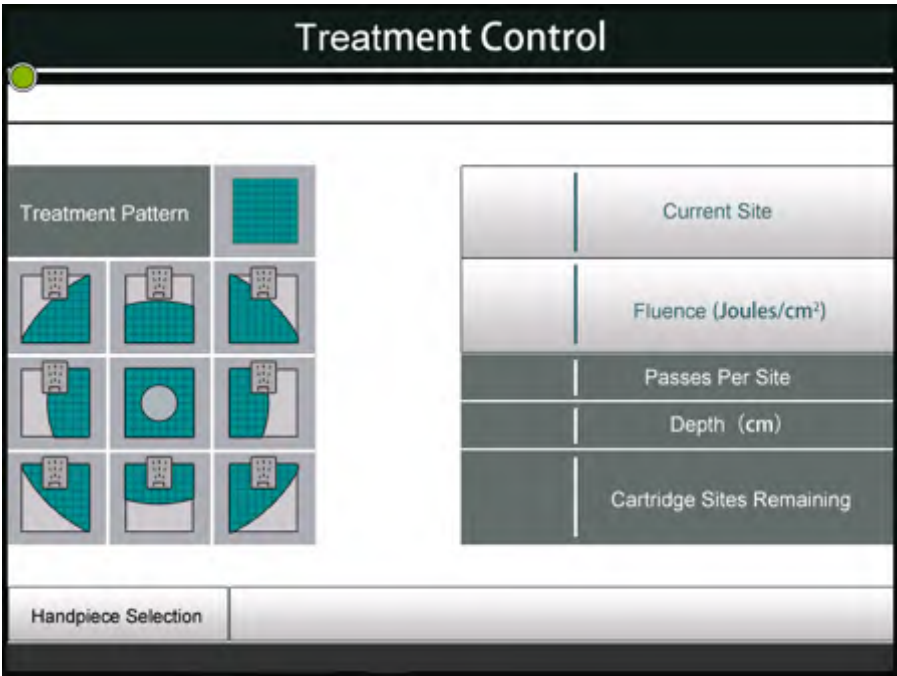

Liposonix



Vaginal

Summary of Advantages

- 

Procedure Time
30 - 45 Minutes
- 

Back To Work
Immediate
- 

Multifunctional
5 handles optional
- 

Full recovery
Immediate
- 

Software Perloaded
All functions' software perloaded already
- 

Upgrade
By inserting new handles only
- 

Number of Treatments
1 every 12 - 18 months
- 

Risks & Complications
Possible mild redness, swelling and tenderness post treatment
- 

Duration Of Results
12 - 18 months



Before and after



Before After



Before After



Before After



Before After



Before After



Before After

Specifications

Power	1200 W
Frequency	4 ~ 7MHz
Net Weight	25kg
Dimension	400*350*1060mm
Voltage AC	100-240v 50/60Hz

Contact us

