microderma pro VRF1

NEWEST TECHNOLOGY ADVANCED

Vacuum Microneedling

Face Lifting

Skin Tightening

Wrinkle Reduction

Pore Reduction

Acne Scars

Scars

Stretch Marks



Treatment Range

FACE: Fine lines, Deep Wrinkles, Appearance of large pores, Acne scars



Fine lines



Deep wrinkles



Large pores



Acne scars



BODY: Stretch Marks, Injury or surgery scar,
Loose and Sagging Skin



Stretch marks




Injury or surgery scar



Loose & Sagging Skin

How It Works

01

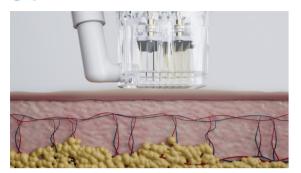


The Micro-needle Cartridge makes perfect contact with the skin

02




The skin in the treatment area is sucked up by the vacuum of the handle

03

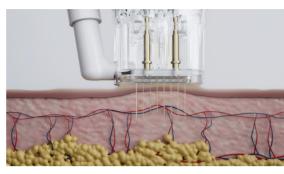


Micro-needles smoothly penetrate the skin with less pain

04

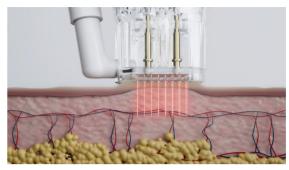

Bipolar RF energy denaturalizes the tissue

05

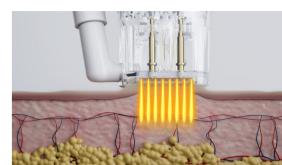

Collagen regeneration and new elastic production begins









Selectable Needle Type

NON-INSULATED NEEDLES

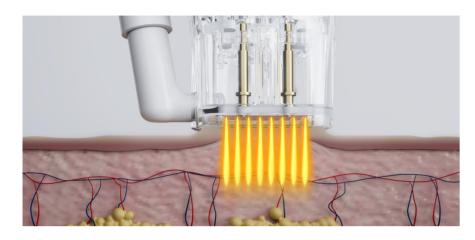


Non-insulated microneedles emit Radio frequency energy over their entire suraceare allowing both the upper and lower layers of the skin to be heated all at once. Non-insulated needles are believed to give more energy due to the increased heating they apply to several layers of the skin at once.

INSULATED NEEDLES

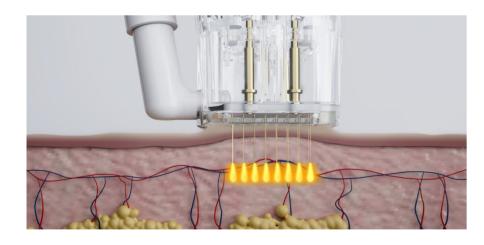

Insulated needles deliver energy only at the tip of the needle, insulated needles will give less energy if compared with non-insulated needles. But Insulation means that the epidermis is spared, no heat energy is dispersed here, only the dermis for a more targeted effect.

Vacuum Rf Microneedling Handpiece

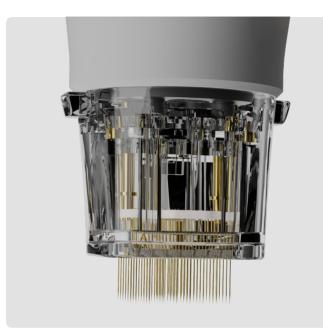

The Microneedling handpiece equipped with USA Haydon stepping motor offers smooth needle delivery that enables comfortable treatment for both practitioners and patients.

MULTIPLE NEEDLE CARTRIDGES

10pins 25pins 36pins 49pins 64pins, needle Depth 0.1mm-4mm adjustable



10 PINS



25 PINS



36 PINS




49 PINS




64 PINS





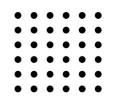

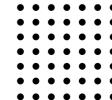

 The identification chip and a sensor inside the handle can automatically detect the needle type and needle output depth,ensuring the safety and accuracy of treatment



USA Haydon stepping motor,
 working life up to 10 years







GOLD PLATED NEEDLES

Surgical steel and plating products have excellent durability, while needles prevent allergic reactions.

0.18MM THIN NEEDLES


The diameter of the needle is only **0.18mm**. Compared with most 0.3mm needles on the market, the thin needle can reduce epidermal and subcutane—ous trauma and speed up skin repair.

The needle tip is only 0.01mm

Thermagical RF Handpiece



The non-invasive Thermagical RF Handpiece is used for skin tightening, anti-wrinkle, lifting, shaping and even fat reduction treatments. No downtime and recovery time.

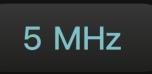

8 MHz

5 MHz and 8 MHz two different frequency output modes, the user can switch the output frequency arbitrarily according to different treatment items and treatment parts, and the RF output depth can be adjusted by switching the output frequency.

FOUR DIFFERENT SIZES


Four different sizes of treatment heads can be applied to the treatment of different parts of the face and the body.



5x5mm



10x10mm

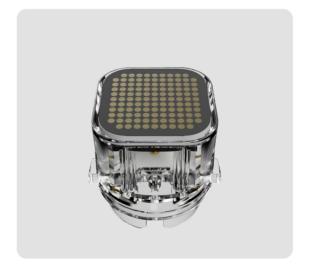

20x20mm




40x40mm


Arms and Legs Are

Rody Area

Cooling Handle

Great for pain reliefing, calming the skin and shrinking pores. Effectively reduce postoperative bleeding and reduce the risk of pigmentation.

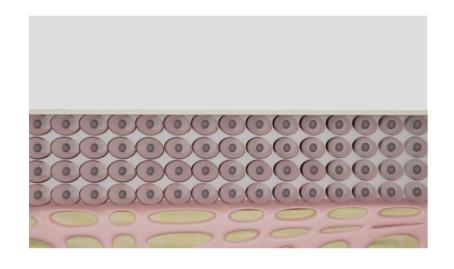

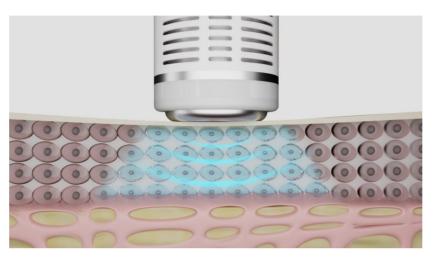





Before And After





Before After



Before After



Before After

Specification

400W Power:

RF Frequency: 5MHz/8MHz

Screen: 15.6 inch Screen:

Pulse Width: Microneedling handle 250-600ms

Thermagical handle 1000-3000ms

Depth: 0-4mm

Net Weight: 22kg

AC220V±10% 50HZ , 110v±10% 60HZ Voltage: